


MED CLAVE

Class N Sterilizer

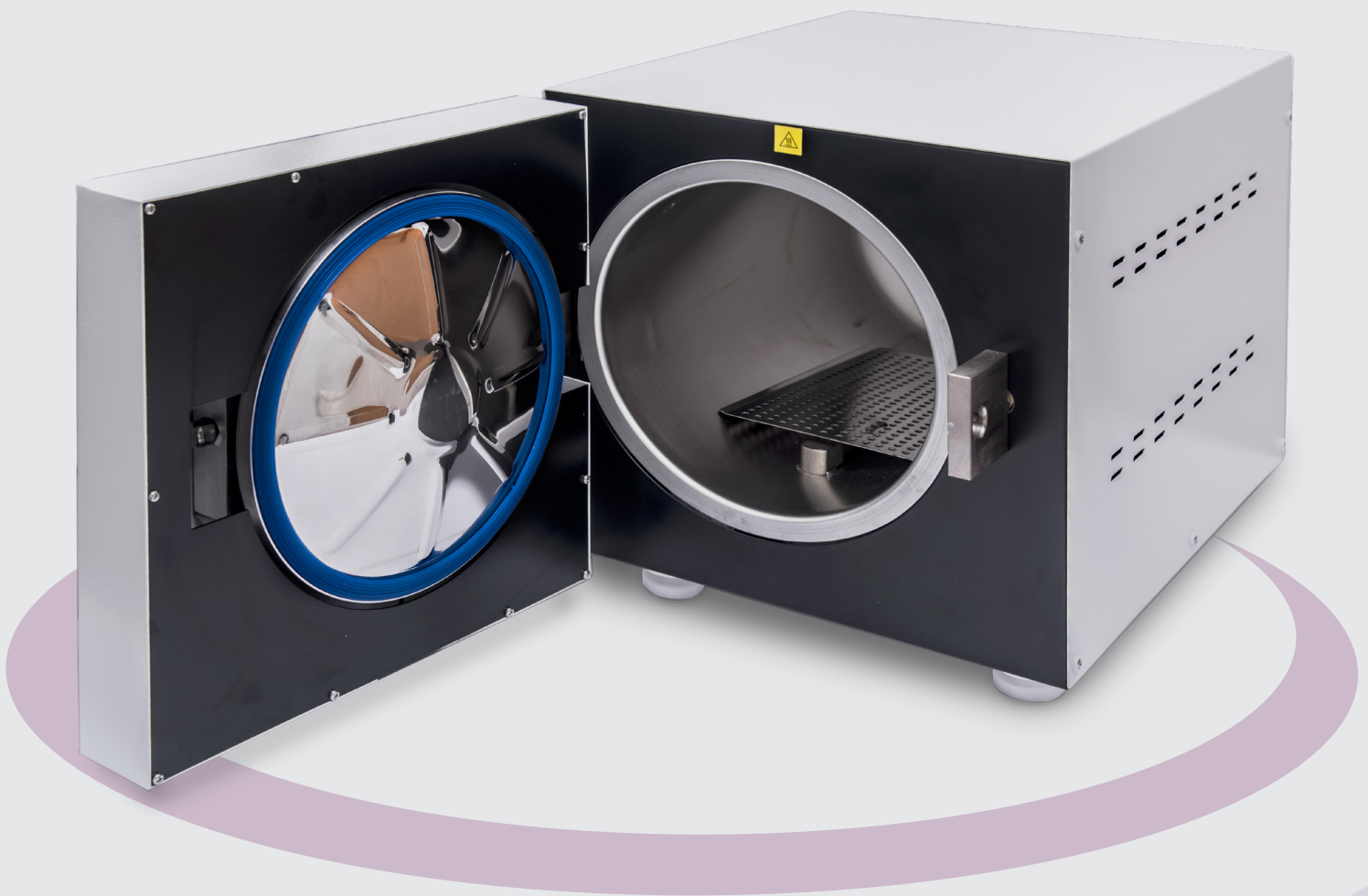

MEDTECH
REMEDIES PVT. LTD.

Safe & Secure

The Med Clave Class N Sterilizer is a modern and versatile product that provides highly effective and safe sterilization for instruments with a 23 Litre capacity. The Med Clave N Class is integrated with microprocessors and sensors for automatic checking and precision control. It has a LED screen panel with a friendly user interface to identify programs and comes with 2 preset sterilization cycles. A desktop model that is very easy to use and maintain and can be operated in the operational area with steam discharge being collected in a waste water tank.



Sterilization of all wrapped and unwrapped products, nonporous, and solid instruments.

Sterilization temperature 121°C and 134°C with
Sterilization time 4-60 mins.

Auto Alarming and malfunction-indication

Safety valve to ensure the chamber can release the
pressure safely in case of emergency





304 AISI steel, without soldering amendments, the finish of which facilitates internal cleaning and maintenance.



Storage tanks for fresh and used water



Double doors locking-system- mechanical lock with electrical lock ensures that the door will open only at 0 pressure.

Technical Specifications

Model	Class N
Volume	23 Liters
Sterilization chamber size	250*450 mm (Dia. X Length)
Trays	03 nos. wire mesh trays
Door lock	Double safety door lock
Gasket	Silicone rubber gasket
Display	LED display
Tanks	4Ltr – 2 Nos
Sterilization time	4 – 60 minutes
Max. pressure	240 Kpa
Max. Temp	137°C
Sterilization temperature	121°C & 134°C adjustable
Product Size	38.5 x 38.6 x 57.6cm
Net Weight	40Kgs
Safety features	Low water & over heat protection, auto-alarm, malfunction indications
Power supply	220 Volts 50/60Hz

For more details:

+91 8590024281

india@med-tech.com

med-tech.in

@medtechremedies